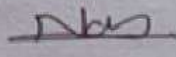



13	Asset Value	2827
14	* Percentage value of the cost of spares with respect to Cost of Purchase/ Asset Value	65.08%
15	Abstract of Service Report provided by the OEM/ Authorized Service Provider/ CYRIX (Attached or Not)	Cyrix
16	Reasons for recommending the equipment as BER	Machine aged upto 5 years. Spare costs are higher. Economically not repairable. So recommending BER.
17	Name & Signature of CYRIX Authority	Nidhish Roy 


*Not mandatory #Based on the period of life and value as per the BER guidelines

* Attach Photograph

Remarks and Recommendations of Junior Consultant (Biomedical) NHM:


Asset value - 2827
 Cost of spares - 1840 (65.08%)
 Life period - 5 years 10 months
 As per the BER guideline recommended as BER

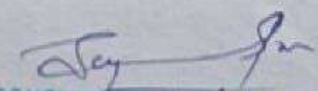
 CONSULTANT BIOMEDICAL ENGINEER
 NATIONAL HEALTH SERVICE
 KOTTAYAM
 Signature of JC BM (NHM)

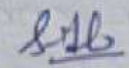SINDHU.K.K 
 Senior Nursing officer.

Date




 Signature of Superintendent / Medical Officer (i/c)

RMO
 General Hospital
 Pala





ELKOMETER
S.N.O. 685215

MIRA
MIRA Calibration Centre
MIRA No. 685215
Cal. on 24.09.11
Exp. on 23.09.13

TANK LOCK OPERATING INSTRUCTIONS
BE OPENED THE APPROPRIATE TANK LOCK HANDLE OF THE LOCK TOWARDS FRONT AND BY ROTATING HANDLE IN THE PLANE OF THE TANK TO CLOSE THE APPROPRIATE TANK LOCK HANDLE



ELKOMETER®
IS 73390


CML-0983369
SPHYGMOMANO
METER 300
MFD. BY
ANITA INDUSTRIES
29, BHAGAT SINGH MARKET
NEW DELHI - 1
S.NO. 315786

MES
Medical Engineering & Services
Thiruvananthapuram - 21. Ph: 0487-2420987
MES CALIBRATION CENTRE
ID: SI No: 68052175
Cal Date: 24.09.21
Exp Date: 23.09.22

TANK LOCK OPERATING INSTRUCTIONS

TO OPERATE THE APPARATUS
TURN THE HANDLE OF THE
LOCK TOWARDS RIGHT
AND THE MERCURY WILL
RISE IN THE GLASS TUBE.
TO CLOSE THE APPARATUS
TILT THE

APPARATUS TOWARDS
LEFT TILL THE COMPLETE
MERCURY GOES IN THE
TANK AND SIMULTANEOUSLY
TURN THE HANDLE OF THE
MERCURY TANK TOWARDS LEFT
NOW BRING THE APPARATUS
BACK TO THE NORMAL POSITION

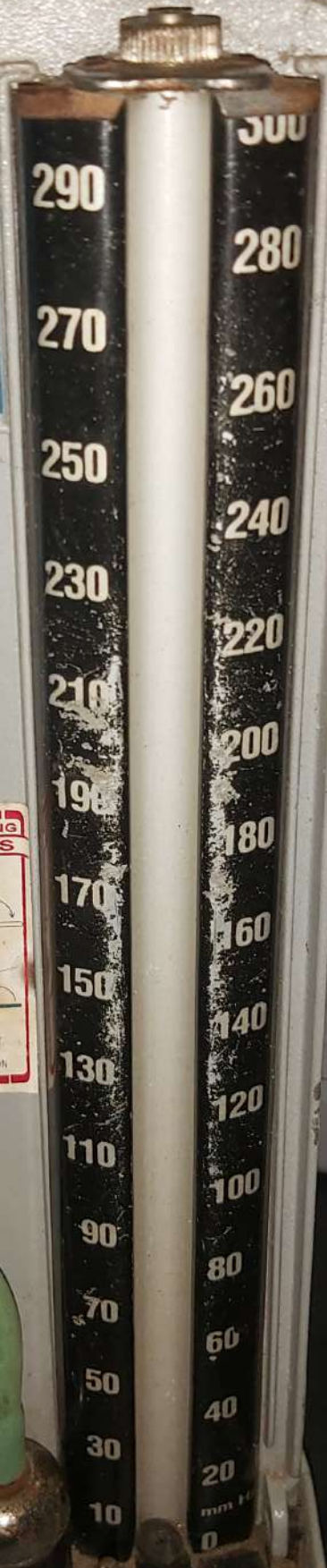

OFF ON

94

BP Apparatus

		Rd	Date of receipt
9/10/18	Rd. from the store.	1	
4/2/19	RECEIVED FROM HIC (L.P)	02	4/2/19
6/5/19	Received from store (DIGITAL BP APPARATUS)	2	
2/3/2020	Rd. from store. (Mercury) (Blue colour) free	1	
	Condensed on per sub. cu		
30/4/22	Tr is from ESTIC (Digital)	1	

SS. Page

Ad-d r

Apparatus

Rd	Date of receipt	cost	RB	78k	Received
1		2804/-	8	location	Bit
02	4/2/19	2332/- costs	10		Amo
2		2909.50	12		kel
1			13		sp
1			3	10	slp 31



**BIOMEDICAL EQUIPMENT
MAINTENANCE PROGRAMME
UNDER
NATIONAL HEALTH MISSION**



SERVICE PROVIDER
Tender No. WO-37/2021-2022/698

CYRIX
HEALTHCARE PVT LTD

No. : 170655

ISO 13485 : 2012 & ISO 9001-2008 CERTIFIED COMPANY | AERB Approved Service Agency

Service Report

30/64 1 B, Petta Junction, Poonithura, Kochi - 682 038, Kerala
Ph : 98472 99500 Website : www.cyrix.com | E-mail : bemp.kl@cyrix.in

Health Facility <u>GH Pala</u>	Call Registration Date : <u>19-12-2003</u>
Address <u>Pala</u>	Caller ID : <u>86514</u>
..... <u>Kottayam</u>	Date of Visit : <u>22-12-2003</u>
Ph : <u>8113968107</u>	Asset No. : <u>0511095</u>
	EQPT Name : <u>BP Apparatus</u>
	Manufacture <u>Anita</u> Model : <u>Ellameter</u>
	S. No. <u>315786</u> Dept. <u>Casualty</u> ..

Service Classification : Breakdown Call PMS Calibration Cust.Training

Problem Identified : Glass tube damaged, Mercury lost, BP cuff nil, Bladder damaged, BP connector required.

Action Taken : Checked. Machine is totally damaged conditions. Required BP cuff, bladder, BP bulb, BP connector, BP glass tube and mercury. Economically not repairable.

Completed Date : 22-12-2003 Time : 01:00pm Spare Required

Spare Replaced Requested

Description	Qty.	Part Number	PR Number
1.			
2.			
3.			

Cyrix Engineer	Date	Start Time	End Time
<u>Nidhish Roy</u>	<u>22-12-2003</u>	<u>12:30pm</u>	<u>01:00pm</u>

Customer Remark Completed Pending
Machine fully in damaged condition.
SINDHU-KU Senior Nursing officer.

22/12/2003
 General Hospital
 Pala

Service Engineer Name : Nidhish Roy
 Signature : [Signature]
 Date : 22-12-2003
 Contact Number : 9072575296

Customer Name :
 Signature :
 Date :
 Contact Number :
 Designation :
 Hospital Seal :

**REPAIR OF BIOMEDICAL EQUIPMENTS UNDER BIOMEDICAL EQUIPMENT
MAINTENANCE PROGRAM (BEMP)**

Recommendations for Beyond Economic Repair (BER)

PROFORMA

Sl. No	Particulars	Details
1	Name of District	Kottayam
2	Name of Hospital	GH Pala
3	Name of Equipment with Make, Model and Serial Number	BP Apparatus Make: Anita Industries Model: Elkometer Sb: 315786
4	Equipment ID & Barcode	0511295
5	Date of purchase / Year of manufacture / Installation Date	09-02-2018
6	Warranty details (Yes/No)	No
7	*AMC/ CAMC Period agreed at the time of purchase	No
8	Date of breakdown (Date of registration of complaint through email/ Toll free)	19-12-2023
9	Action taken	Checked & serviced
10	Present status of the equipment (Fully damaged / partially damaged)	Fully damaged.
11	Recommendations for repair (required service details)	Glass tube damaged. Mercury lost BP cuff nil, BP glass tube required, BP bladder damaged BP connector required.
12	Cost of spares (specify parts and cost)	Glass tube - 200 Mercury - 1200 BP bladder - 180 BP bulb - 60 BP connector - 20 BP cuff - 180, Total = 1840